



## **Srl<sub>2</sub> Scintillation** Packaging



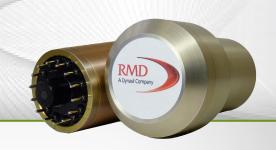

1.5" x 1.5" Srl, Crystal

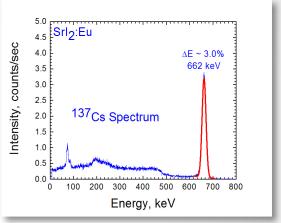

1.5" Srl, Crystal - Energy Spectrum

Srl<sub>2</sub>

RMD supplies Srl, scintillation crystals designed to meet our customer's needs. Srl, scintillation crystals and packaging are specified after discussions with RMD and are quoted based upon the customer's specification. Available packaging includes:

**Unpackaged: Shipped in an inert atmosphere container.** Srl<sub>2</sub> scintillation crystals are hygroscopic. Unpackaged Srl<sub>2</sub> crystals are available for customers who have the ability to work with Srl<sub>2</sub> in a dry box and who are proficient in crystal preparation and packaging.

Hermetic Package: With an optical window. Srl<sub>2</sub> scintillation crystals are supplied in a right-circular cylindrical package with an optical window on one end. The housing protects the right cylindrical shaped crystal that is permanently sealed inside the housing.

**Packaged with a Photomultiplier Tube.** Srl<sub>2</sub> scintillation crystals are packaged in conjunction with a suitable photomultiplier tube, or other photodetector.

Other crystal shapes, sizes and packages are available upon request.

Model Number	Scintillation Material	Crystal Shape and Crystal Dimensions	Energy Resolution	Package Style
SRI-25-PHI-25-S-117	Srl <sub>2</sub>	Cylindrical 25 mm diameter, 25 mm length	< 4 % FWHM at 662 keV line of <sup>137</sup> Cs	Sealed Package
SRI-25-PHI-25-P-121	Srl <sub>2</sub>	Cylindrical 25 mm diameter, 25 mm length	< 4 % FWHM at 662 keV line of <sup>137</sup> Cs	Scintillation Crystal is permanently mounted to a photomultiplier tube.
SRI-38-PHI-38-S-117	Srl <sub>2</sub>	Cylindrical 38 mm diameter, 38 mm length	< 4 % FWHM at 662 keV line of <sup>137</sup> Cs	Sealed Package
SRI-38-PHI-38-P-121	Srl <sub>2</sub>	Cylindrical 38 mm diameter, 38 mm length	< 4 % FWHM at 662 keV line of <sup>137</sup> Cs	Scintillation Crystal is permanently mounted to a photomultiplier tube.

